GERMAN
POOL

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Electronic Body Scale
人體電子磅

ESB-115

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Please read these instructions and warranty information carefully before use and keep them handy for future reference.

使用說明書 USER MANUAL

目錄 Table of Contents

注意事項 Notice Before Use 03

操作說明 Operating Instructions 04

清潔與保養 Care & Maintenance 05

常見故障及處理 Troubleshooting 05


技術規格 Technical Specifications 06

保用條款 Warranty Terms & Conditions 07

注意事項 Notice Before Use


1. 為了使您更好的使用本產品，我們建議您認真閱讀本說明書，並妥善保存以供參閱。
 2. 本產品只供家庭測量體重，不可用於商業或其他用途。
1. In order to ensure you have a better understanding of the product, we recommend that you read the manual carefully and retain it for reference.
 2. This scale is only used in household for body weight measurement, and not for commercial or other use.

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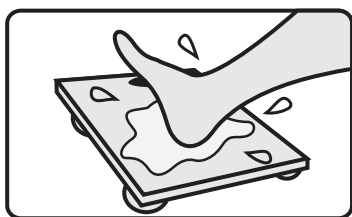
- 請勿濕腳磅重，謹防滑倒。
 - 當玻璃平台被弄濕後會變得很滑，所以必需要保持玻璃平台乾燥。
- Do not weigh with wet feet. Beware of slipping.
 - The glass platform will be slippery when it is wet. Keep the glass platform dry.

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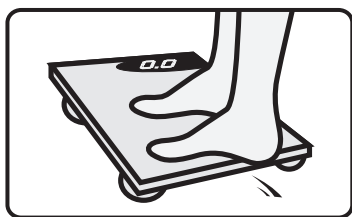
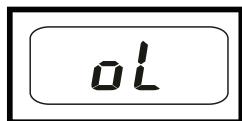
Overload Indication: Scale overloaded. Step off from the scale to avoid damages to the product itself. (Max weight capacity is indicated on the scale).

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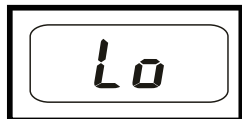Low Battery Indication: If the voltage of battery is lower than $2.4V \pm 0.2V$, the screen will display "Lo". Replace new batteries.

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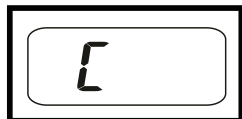
Re-Test: Displaying "C" when you step off the scale means that an error has occurred. Re-measure your weight again.

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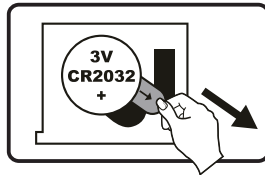1. Removing the plastic tab from the battery before using the scale

- The scale includes one 3V lithium battery (CR2032).
- Remove the screw and open the battery compartment on the back of the scale. Gently pull out the plastic tab to activate the battery.
- Close the battery compartment.

Warning: If you are not using the product for a long time, we recommend that you remove the battery.

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
2. Placing the scale


Put the scale on a flat and rigid ground before use.

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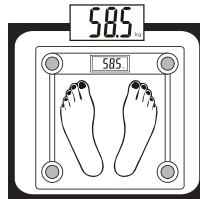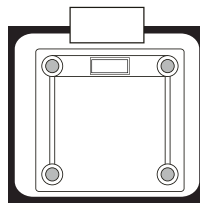3. Using the scale

- Step on the scale. Please remain stationary until your weight shown on the display is stable.
- The screen displays your weight value. It flashes for 3 times and is then locked. The value is then your body weight.
- If you want to measure again, step down from the scale and then step back on it.
- Scale will automatically shut down after no use for a while.

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- If the scale fails to power on

Please check whether the battery is installed properly or it has run out of power. Insert or replace a new battery if necessary.

技術規格 Technical Specifications

型號 Model	ESB-115
高精準的STRAIN GAUGE稱重技術 High-precision STRAIN GAUGE weighing technology	有/YES
節能液晶顯示屏幕 Energy-saving LCD screen	有/YES
自動開關機功能 Auto on/off function	有/YES
自動零點重設 Auto zero resetting	有/YES
低電量顯示 Low Power Indication	有/YES
超載指示 Overload Indication	有/YES
最大磅重量 Maximum weighing capacity	150kg / 330lb
磅重量範圍 Loading range	5kg (11lb) ~ 150kg (330lb)
稱量刻度 Weight Graduation	100g , 0.21b
誤差 Weight Deviation	±0.75%+ID
適用環境 (溫度) Operation temperature	10-40 °C
淨重 Net Weight	1.5kg
電源 Battery	1枚3V CR2032電池 (附送) 1 x 3V lithium battery (CR2032) (included)

產品規格及設計如有變更，恕不另行通知。
Specifications are subject to change without prior notice.

中英文版本如有出入，一概以中文版為準。
If there is any inconsistency or ambiguity between the English version and the Chinese version,
the Chinese version shall prevail.

本產品使用說明書以 www.germanpool.com 網上版為最新版本。
Refer to www.germanpool.com for the most up-to-date version of the Operating Instructions.

請填妥以下空格並與發票一併保存，以便維修時核對資料。用戶提供本公司之所有資料必須正確無誤，否則免費保用將會無效。
Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty and the blade therein a lifetime warranty provided by German Pool (effective starting from the date of purchase).

* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online:
www.germanpool.com/warranty.
- 2) Complete all fields on the attached Warranty Card and fax it back to our Customer Service Centre along with the purchase invoice.

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
 - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
 - damages caused by accidents of any kind (including transport and others);
 - blender jar, plastic parts and accessories;
 - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
 - purchase invoice is modified by unauthorized party;
 - product is used for any commercial or industrial applications;
 - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
 - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

Customer Service & Repair Centre

German Pool (Hong Kong) Limited
8/F, China Travel (Cargo) Logistics Centre
1 Cheong Tung Road
Hung Hom, Kowloon, Hong Kong
Tel : +852 2333 6249
Fax : +852 2356 9798
Email : repairs@germanpool.com

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This warranty is valid only in Hong Kong

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Room 113, Newport Centre Phase II, 116 Ma Tau

Kok Road, Tokwawan, Kowloon, Hong Kong.

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© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

www.germanpool.com