GERMAN
POOL

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Far Infrared Tapping Massager
遠紅外線肩部按摩器

MS100

在使用之前請詳細閱讀「使用說明書」及「保用條例」，並妥為保存。
Please read these instructions and warranty information carefully
before use and keep them handy for future reference.

使用說明書 OPERATING INSTRUCTIONS

Table of Contents

Notice Before Use	03
Product Structure	05
Operating Instructions	07
- Wearing Method	
- Operation	
- Mode Selections	
Care & Maintenance	10
Technical Specifications	10
Warranty Terms & Conditions	19

Notice Before Use

Important

For safety purpose, you should read the instructions carefully before using the appliance for the first time. Save the instructions for future reference.

The appliance is designed for domestic use and may only be installed and operated in accordance with this instruction manual.

1. Do consult your doctor before use if you have any concerns about your health.
2. Pregnant women, or people suffering from heart disease, hypertension, or vertebrate illness, should use this product under the guidance of a doctor.
3. Mobility-handicapped persons are not recommended to use this product.
4. Close supervision is necessary when this appliance in used by, on, or near children, invalids, or disabled persons.
5. Do not place this product where children are close.
6. This product is not a professional medical equipment, not specifically designed to relax the muscles, not for any medical purposes.
7. This product should not be used on any swelling part.
8. Do not use this product on sensitive skin, and when inflammation happens, stop using this.
9. Do not apply Warp Neck & Shoulder Massager directly on head, joint of elbows and knees, breasts or abdomen.
10. Do not cover this product with any articles. Do not cover the venting holes.
11. Do not use this product on bed and never operate the appliance with the air opening blocked, such as under blankets, pillows, which may result in overheating and cause fire, electric shock or injury to users.
12. Do not use this products when in the bath or under water. This product is not for use in damp environment.
13. If this product falls in water by accident, please turn off the power and unplug the power cord immediately. Taking it out directly from water is prohibited.
14. Do not expose the product to heaters, stoves or other source of heat.
15. You are advised to use this product 1 hour after meals, and remember to drink 150 cc water beforehand.
16. Another 250 cc of drinking water is recommended after massaging.
17. You are advised to use this product for less than 30 minutes each time.
18. Do not press this product hard with foot.

Notice Before Use

19. Long time of using this product may cause overheat, which may shorten the life of the product. If this happens, stop using it immediately.
20. You are not advised to use this product before going to bed, because this product may boost your spirit.
21. Unplug the power cord when the product is not in use for a long time.
22. Before cleaning this product and after use, turn off the massager and unplug the power cord by holding the plug, not the cord.
23. Do not use this product in the outdoors and in the room where temperature is over 40°C.
24. To turn off this product, turn all controls to the off position, then remove plug from outlet.
25. This product is developed for personal use and not for commercial purposes.
26. Do not wind the power cord around the product. Do not operate the product with a damaged cord or tampered cord, as this may damage the cord or cause electric shock.


Product Structure

Venting Holes
(Caution: Do not block venting holes.)



Product Structure

Control Panel

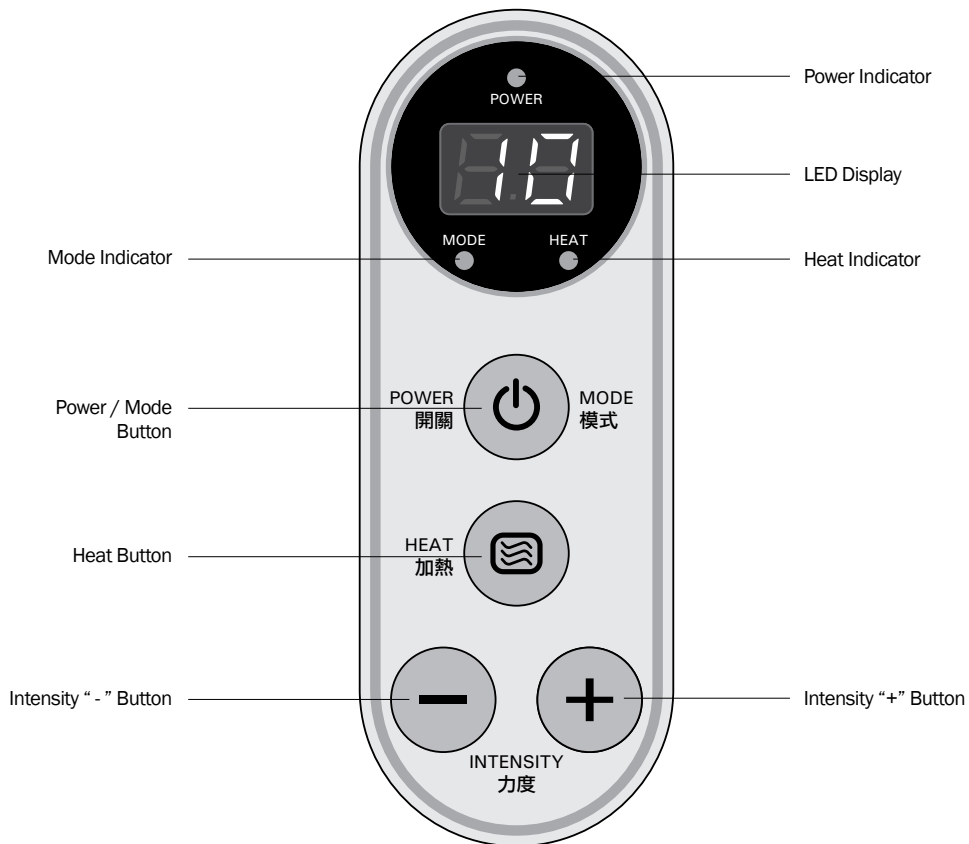Operating Instructions

A) Wearing Method

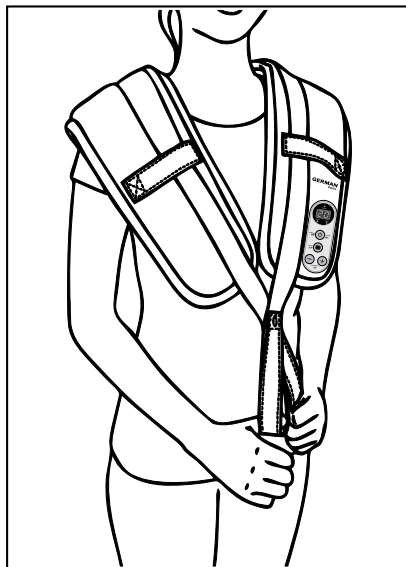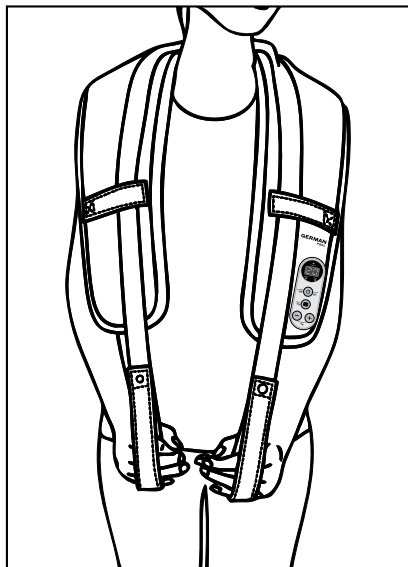Control panel facing outward, apply the product over the body part to be massaged. You may hold the handle strap to secure the position or fasten the button on the handle.


Neck and shoulder


Back

Lower back

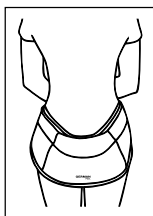Waist

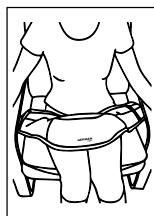Thighs

Operating Instructions

B) Operation

Note:

Connect the power plug, and turn on the main switch, place the massager on the body part and make sure it is comfortably secured.

Automatic Mode (Blue Light)

1. Connect the power. A long beep is heard, and the Power Indicator will flash in red. Press Power/Mode button  for about 2 seconds, the massager will enter Automatic Mode; The Mode Indicator will light up in blue.
2. Massager will start its operation from the first massage mode. It automatically switches to the next mode until the 20th mode and cycles again (Please refer to "Mode Selections" on next page). The default operation time is 8 minutes.
3. a. Setting Heating Function: Pressing the Heat button  will start the heating function. Heat Indicator lights up.
b. Adjusting Massaging Intensity: The default intensity is 1. Press Intensity **+** or **-** button to increase the intensity level, (each press will be accompanied by a short beep), and the display will show the selected intensity level.
4. Press and hold Power/Mode button  for 3 seconds to turn the power off.

Manual Mode (Red Light)

1. Under Automatic Mode, press Power/Mode button  to switch from Automatic Mode to Manual Mode. The buzzer will give a short beep and the Mode Indicator will turn from blue to red.
2. Each press of Power/Mode button  will switch the massager to the next massage mode (Please refer to "Mode Selections" on next page).
3. a. Setting Heating Function: Pressing the Heat button will start the heating function with indicator light on.
b. Adjusting Massage Intensity: Press Intensity **+** or **-** button to adjust the intensity level, each press will be accompanied by a short beep, and the display will show the selected intensity level.
4. Press Power/Mode button  for 3 seconds to turn the power off.

Operating Instructions

C) Mode Selections

Mode	Massage Actions
1	Switch tapping
2	Resonance tapping
3	Fast switch tapping
4	High frequency switch tapping
5	Variable-frequency switch tapping
6	High-frequency double vibration
7	High-frequency crisscross tapping
8	Intermittent tapping
9	Variable-frequency double vibration tapping
10	Alternating tapping
11	Fast surrounding tapping
12	Surrounding crisscross tapping
13	Mild crisscross tapping
14	Accelerated crisscross tapping
15	Continuous crisscross tapping
16	Left-right rapid crisscross tapping
17	Right-left rapid crisscross tapping
18	High-speed double vibration tapping
19	Left-right mild tapping
20	Right-left soft tapping

Note:

The product is equipped with an overheat protection device. Long hours of continuous usage may cause the product to stop operation due to overheating. When this happens, allow the product to cool down before using again.

Care & Maintenance

- Before cleaning this product and after use, please power off the massager and unplug the power cord by holding the plug, not the cord.
- Do not use this products when in the bath or under water. This product is not for use in damp environment.
- If this product falls in water by accident, please turn off the power and unplug the power cord at once. Do not remove it from water immediately.

Technical Specifications

Model No	MS100
Voltage / Frequency	220 V ~ 50 Hz
Rated Power	100 W

Specifications are subject to change without prior notice.

If there is any inconsistency or ambiguity between the English version and the Chinese version, the English version shall prevail.

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Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online: www.germanpool.com/warranty.
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
 - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
 - damages caused by accidents of any kind (including transport and others);
 - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
 - purchase invoice is modified by unauthorized party;
 - product is used for any commercial or industrial applications;
 - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
 - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

Repair Centre - Small Appliances

German Pool (Hong Kong) Limited
8/F, China Travel (Cargo) Logistics Centre
1 Cheong Tung Road
Hung Hom, Kowloon, Hong Kong
Tel : +852 2215 4582
Fax : +852 2355 7100
Email : repairs@germanpool.com

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This warranty card is valid only in Hong Kong

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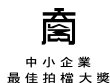
© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

www.germanpool.com