GERMAN
POOL

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eye care massager | me10

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在使用之前請詳細閱讀「使用說明書」及「保用條例」，並妥為保存。
Please read these instructions and warranty information carefully before use and keep them handy for future reference.

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IMPORTANT

For safety purpose, you should read the instructions carefully before using the appliance for the first time. Save the instructions for future reference.

The appliance is designed for domestic use and may only be operated in accordance with these instructions for use.

- Persons who have had eye surgery or suffer from glaucoma, cataract, retinoschisis or other serious eye diseases should not use the massager.
- Take off contact lenses, hairpins and hair accessories before using the massager.
- Do not wash the device with water. If cleaning is needed, wipe it with aseptic cotton dipped with alcohol.
- Hold the main unit and do not only hold the cable of the device.
- Do not disassemble the device by yourself.
- Since the cable is on the left side of the massager, you may feel slightly heavier on one side. This is a normal phenomenon and does not affect its function.
- Put on the eye massager before turning it on. This can prevent damaging the air bags.
- Do not use the device in an environment with high temperature or humidity.
- Avoid strong sunlight and corrosive gas.

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The eye massager functions through vibrating, kneading and heating. Adopting state-of-the-art microcomputer system and far infrared technology, it can improve blood circulation of eyes, help increase oxygen supply, adjust the function of ciliary muscle, relax the eyes and brain, promote metabolism and maintain healthy eyesight. The device also adopts different relaxing music, α brain wave music and soothing natural sounds to stimulate brain cell activities.

Scope of Application


It is mainly used to relax your eyes and brain. Meanwhile, the massager can help to prevent myopia, pseudomyopia, lengthening and distortion of optic axis.

Target Users

- People who suffer from myopia, amblyopia, astigmatism, such as students.
- People who use the eyes frequently during work, including computer operator, accountants, designers and engineers.
- People who suffer from pouch, dark eye circle and eye wrinkles caused by tiredness and clogged blood circulation.
- People suffering from insomnia, headache and presbyopia.

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
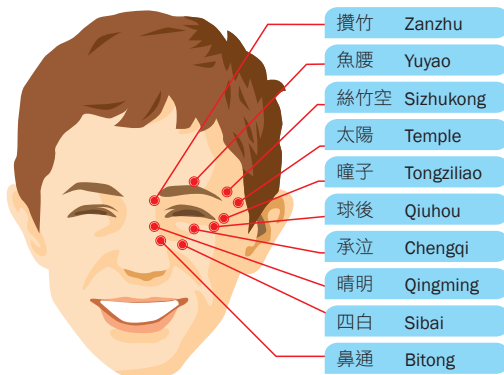
Zanzhu	Headache, eye gall, blurred vision, walleys
Yuyao	Keratitis, facial paralysis, eye ophthalmoplegia
Sizhukong	Hemicrania, remove obstruction in channels to clear eyesight
Temple	Headache, eye disease
Tongziliao	Refractive error, optic atrophy
Qiuhou	Dispel wind and remove heat, moisten and clear eyes
Chengqi	Acute and chronic keratitis, epiphora induced by wind, retinitis
Qingming	Dispel wind and remove heat, moisten and clear eyessuralgia
Sibai	Facial paralysis, trifacial neuralgia
Bitong	Refractive error, nastitis and nasal congestion

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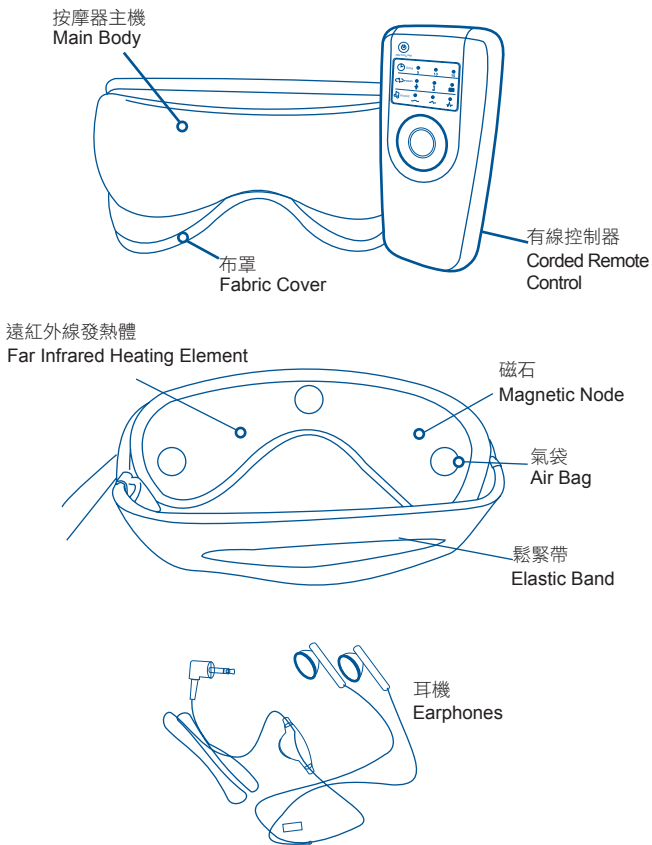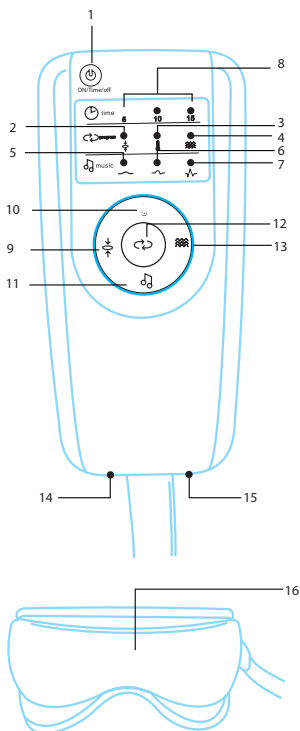1. ON/OFF/TIMER
(Press to turn on, press again to select time, press and hold to turn off)
2. Air Kneading Indicator
3. Heating Indicator
4. Vibrating Indicator
5. Relaxing Music Indicator
6. α Brain Wave Music Indicator
7. Natural Music Indicator
8. Timing Indicator (15 min/10 min/5 min)
9. Air Kneading Mode Selection
10. Heating Mode Selection
11. Music Mode Selection/Switch
12. Massaging Mode Selection
13. Vibrating Mode Selection
14. Power Jack (DC 6V, 1000 mA)
15. Earphone Jack
16. Working Indicator
(Light flashes during operation, and stops flashing when operation stops)

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
Battery Installation

1. Slide up and open the battery compartment cover at the back of remote control.
2. Insert four AA batteries according to the polarity marks inside (Alkaline batteries are recommended).
3. Put on the battery cover.

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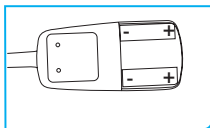Connecting Power Adapter (Purchased Separately)

1. Use DC6V adapter.
2. Insert the adapter to the power jack on the remote control.
3. Connect the adapter to power.

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
1. 將耳機插入機身耳機插孔。
Connect earphones to the socket.

2. 請將布罩對準雙眼戴上主機。
Position and align device with eyes.

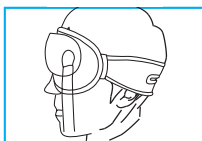3. 調整鬆緊帶扣至舒服位置。
Adjust the elastic band to comfortable length.

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Start Up Status & Default Setting

- Press On/Off button, the device starts to work.

- Working Indicator flashes
- Air Kneading Indicator on
- Heating Indicator on
- Vibration Indicator on
- 5-min Timer Indicator on
- Music: Relaxing music
- Air Kneading: massage temples and eyes

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Select Functions

1. Music

- Press Music button to choose music or no music
- Press Music button shortly, you can choose any music mode.
 - Relaxing music
 - Brain wave music 1
 - Brain wave music 2
 - Natural music (dripping)
 - Natural music (bird singing)
 - Natural music (ocean wave)
 - Natural music (cricket)
 - Natural music (flowing water)
- Press and hold Music button to switch off the music function.

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2. Air Kneading

- Press Air Kneading button, massage setting will change in the following order:
 - Massage Temples + Eyes
 - Massage Eyes
 - Massage Temples

3. Heating

- Press Heating button, heating setting will change in the following order:
 - Low Temperature
 - High Temperature

4. Vibrating

- Press Vibrating button, vibrating setting will change in the following order:
 - Strong/Medium/Weak/ Auto
 - Comfortable/Press (Heavy)/ Press (Medium)/Press (Quick)
 - Alternate/Alternate (Quick)

5. Timer

- Press Timer button, time setting will change in the following order:
 - 5 Minutes
 - 10 Minutes
 - 15 Minutes

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Power Supply Inspection

- Inspect the power supply. Do not turn on the device if the power voltage exceeds the specified range.
- Inspect device and cable regularly. If any damage or breaking occurs, stop using the massager immediately and send it to our Repair Centre.

Maintenance

- Turn off the massager before cleaning the remote control.
- Keep the device at safe places and do not put any heavy objects on top of it.
- Take measures to prevent dampness, mice and pollution.
- Do not disassemble the device by yourself in case of breakdown. Please contact our Repair Centre for repair.

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When you run into any problem operating this product, use the table below to troubleshoot the problem. If the problem persists, contact our Repair Centre.

故障現象 Phenomenon	原因分析 Possible Cause	處理方法 Solution
無法開啟電源 Cannot turn on the device.	電纜是否完好無損？ Check whether the cable is in good condition.	檢查電纜連線狀況 Check the connection of the cable.
	電池電量是否不足？ Check whether the power level of batteries is too low.	如果是，請更換新電池。 If yes, change the batteries.
	電池是否安裝正確？ Check whether the batteries are placed properly.	請確保正確安裝電池電極。 Make sure the polarities of the batteries are correct.
請確保將電池裝入並固定在電池盒中。 Make sure the batteries are installed and fixed in the case.		
指示燈亮度不夠 The indicators are not bright enough.	電池電量是否微弱？ Check whether the power level of batteries is too low.	如果是，請更換新電池。 If yes, change the batteries.
氣壓、加熱強度太弱 The air pressure is low and vibration is weak.	電池電量是否微弱？ Check whether the power level of batteries is too low.	如果是，請更換新電池 If yes, change the batteries.
		請確保將電池裝入並固定在電池盒中。 Make sure the batteries are installed and fixed in the case.

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* Above table enlists the most common problems regarding the operations of the product and their possible causes as well as suggested solutions. It is strictly for users' reference only.

* Should there be any product malfunction or breakdown, DO NOT repair it yourself. Contact our Repair Centre or consult a licensed technician immediately.

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產品規格及設計如有變更，恕不另行通知，並以產品實物為準。
Product specifications and appearance are subject to change without notice.

中英文版本如有出入，一概以中文版為準。
If there is any inconsistency or ambiguity between the English version and the Chinese version, the Chinese version shall prevail.

請填寫以下空格並與發票一併保存，以便維修時核對資料。用戶提供本公司之所有資料必須正確無誤，否則免費保用將會無效。
Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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Warranty Terms & Conditions

Our Repairs & Customer Service Centre will provide rectification of any defect of manufacturing materials free of charge given that such defect is reported within 12 months from date of sale.

* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online:
www.germanpool.com/warranty
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer)

- Customer should always follow the operating instructions. This warranty does not apply to: damages caused by incorrect voltage, improper operations, unauthorized installation or repairs as well as accidents.

- Customer who fails to present original purchase invoice will be disqualified for free repair services.

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Repair Centre

German Pool (Hong Kong) Limited
2/F, Kai It Industrial Building
15 Mok Cheong Street
Tokwawan, Kowloon, Hong Kong
Tel : +852 2215 4582
Fax : +852 2356 9798
Email : repairs@germanpool.com

此保用只適用於香港 | This warranty is valid only in Hong Kong

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