GERMAN  
POOL

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mini 迷你手提按摩器  
massager | mm5

使用說明書 OPERATING INSTRUCTIONS

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---

Please register your warranty information now!

For Warranty Terms & Conditions,

please refer to the last page of these operating instructions.

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Please read these instructions and warranty information carefully before use and keep them handy for future reference.

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Incorporating quintessences of traditional Chinese medicine and adopting vibration and massage functions, this massager keeps meridians and collaterals clear, activates Qi, and improves blood circulation. It helps to ease chronic pain, acute pain, muscle ache and palsy, and can promote blood circulation, exudation of pain-restricting substances and immediate interdiction of pain transmission.

## Scope of Application

Applicable to the whole body, including arms, thighs, waist and abdomen.

## Target Users



Students, computer users, accountants, designers, clerks, teachers and long-hour working personnel at all ages.

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This device is designed according to the curves of human body. It is equipped with 3 massage heads and a high performance motor, achieving the best massage force.

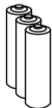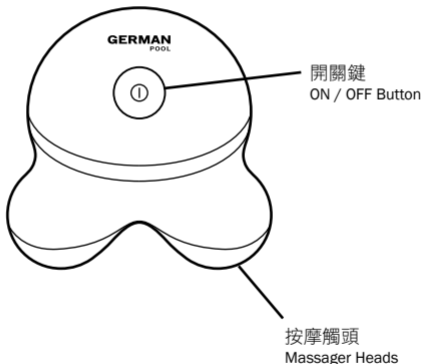隨機附設：乾電池 (AAA x3)  
Items included: Batteries (AAA X3)

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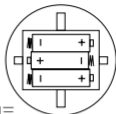
## Battery Placement

1. Open the unit cover in anti-clockwise direction (Fig 1&2). Place three pieces of AAA dry batteries in.
2. Place the battery cover back.

Note: Please check if poles of the batteries are correct (Fig 3) and if battery cover is securely closed.

## Operating Method

1. Place the three kneading massage heads of the massager close to the body part to be kneaded, press ON/OFF button. The machine starts working and the blue LED of the massage heads will light up.
2. Take the main unit away from the body and press ON/Off button, the machine stops working.

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## Applicable Conditions

The massager must be used under the following conditions:

Temperature: 5°C - 40°C

Relative humidity: Max. 80%, without corrosive gas and with good ventilation

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## Portable Usage

The product is equipped with simple vibration-proof measures and suitable for transportation by air, railway and ship. Keep it away from rain, snow, splashing and collision. Keep it upright.

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## Storage

- If storage period is over 6 months, remove batteries from the unit.
- Do not store it close to ground, wall and roof.
- Good ventilation is required. Keep it away from sunshine and corrosive gas.

## 清潔與保養 Care & Maintenance

---

- 清潔前請先關機。
  - 如本機變髒，可用軟布沾上稀釋酒精擦拭機身。
  - 使用後應將按摩器放置於安全地方，切勿在按摩器上放置重物。
  - 採取防潮、防鼠、防污染措施。
  - 當出現故障時應找專業技術人員修理，切勿隨意拆裝按摩器。
- Turn off the unit before cleaning.
  - If the unit body gets dirty, clean it with a soft cloth dipped with diluted rubbing alcohol.
  - Keep it in a safe place after use. It is strictly prohibited to place heavy articles on top of the device.
  - Measures to prevent dampness, mice and pollution are required.
  - In case of failure, please have professional technicians repair it. Do not disassemble it by yourself.

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When you run into any problem operating this product, use the table below to troubleshoot the problem. If the problem persists, contact our Repair Centre.

| 故障現象<br>Phenomenon                     | 原因分析<br>Possible Cause   | 處理方法<br>Solution  |
|--|--|---|
| 按摩器不啟動。<br>The massager does not work. | 電池是否已消耗盡？<br>電池的正負極是否正確？<br>Batteries using out? Correct poles of the batteries? | 更換電池，確認正負極方向，正確插入。<br>Replace the batteries with correct poles. |
| 按摩力度弱。<br>Weak massage force.          | 電池是否消耗盡了？<br>Batteries using out?  | 更換電池。<br>Replace the batteries with correct poles.              |

- \* 以上為常見故障的判斷與檢查，只供用戶作參考之用。
- \* 如遇故障切勿自行修理，避免發生危險及對機件造成損害，請立即找專業技術人員或聯絡本公司作出適當維修。
- \* Above table enlists the most common problems regarding the operations of the product and their possible causes as well as suggested solutions. It is strictly for users' reference only.
- \* Should there be any product malfunction or breakdown, DO NOT repair it yourself. Contact our Repair Centre or consult a licensed technician immediately.

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產品規格及設計如有變更，恕不另行通知。  
Specifications are subject to change without prior notice.

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If there is any inconsistency or ambiguity between the English version and the Chinese version, the Chinese version shall prevail.

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**For additional product information,  
please visit our website.**

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Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

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\* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

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- 1) Visit our website and register online:  
[www.germanpool.com/warranty](http://www.germanpool.com/warranty).
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

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1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
  - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
  - damages caused by accidents of any kind (including transport and others);
  - plastic parts and accessories;
  - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
  - purchase invoice is modified by unauthorized party;
  - product is used for any commercial or industrial applications;
  - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
  - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

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Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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電郵：repairs@germanpool.com

## Repair Centre

German Pool (Hong Kong) Limited  
8/F, China Travel (Cargo) Logistics Centre,  
1 Cheong Tung Road,

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