GERMAN
POOL

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Far Infrared Tapping Massager
遠紅外線肩部按摩器

MS101

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Please read these instructions and warranty information carefully before use and keep them handy for future reference.

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Please register your warranty information now !

For Warranty Terms & Conditions,

please refer to the last page of these operating instructions.

Contents

Warning	04
Product Structure	05 -06
Operating Instructions	07 - 08
Care & Maintenance	09
Technical Specifications	10
Warranty Terms & Conditions	19

Warning

Important

For safety purposes, you should read these instructions carefully before using the massager. Retain these instructions for future reference.

The appliance is designed for domestic use only, and may only be installed and operated in accordance with this User Manual. Please observe basic safety precautions when you operate this appliance.

1. Consult your doctor or a medical professional before using this massager if you have any concerns about your health.
2. Pregnant women, or people suffering conditions such as heart disease, hypertension, and/or vertebrate illness, should only use this massager under the guidance of a doctor.
3. This massager is not recommended for use by the mobility-handicapped.
4. Close supervision is necessary when this massager is used by, on, or near children, invalids, or disabled persons.
5. Do not place this product in an area easily accessible by children.
6. If you feel any physical discomfort during the use of this appliance, or when the appliance fails to function normally, stop using immediately.
7. This product cannot be used for diagnostic or therapeutic purposes, as it is not a professional medical equipment, nor is it professionally designed to relax the muscles.
8. This product should not be used on any swollen part of the body.
9. Do not use this product on sensitive skin. If inflammation results from using the appliance, cease using immediately.
10. Do not apply the Far Infrared Tapping Massager directly on the head, the joints of elbows and knees, or on the breasts or abdomen.
11. Do not cover this appliance. Do not block or cover the venting holes.
12. Do not use this product in bed, or block the venting holes with items such as blankets or pillows during operation. This may result in overheating and risk of fire, electric shock, or injury.
13. Do not use this product when in the bath or under water. The appliance is not meant for use in a damp environment.
14. If you accidentally drop this product into water, DO NOT retrieve it before turning off the power and unplugging the cord.
15. Do not use the product near heaters, stoves, or other sources of heat, especially open flames.
16. It is advisable to drink 0.15L of water before, and 0.5L after using this product. Do not use this product within one hour before or after eating.
17. Do not use this product over a 30-minute continual period. When selecting massage intensity, please adjust according to your personal physical conditions and needs.
18. Do not sit on, or violently stamp on the appliance.
19. Prolonged usage may cause the massager to overheat. This drastically reduces the lifespan of the product. If the product gets too hot, turn it off immediately and wait for it to cool down before using.
20. It is not advisable to use this product before sleep, because the FIR massage may make you more alert.
21. Unplug the power cord after each use.
22. After use and before cleaning, turn off the massager and unplug the power cord via the plug - not the cord. Pulling on the power cord may cause damage to the product.
23. Do not use this appliance outdoors or in rooms with ambient temperature over 40°C.
24. To turn off this product, turn off all the controls before removing the plug from outlet.
25. This product is designed for home use, and is not suitable for commercial purposes.
26. To avoid the risk of short-circuiting or damage to the power cord, do not wind the power cord around the product. Do not operate the product with a damaged cord.

Product Structure


1. Main Unit

*Product specifications are subject to change without prior notice. Please refer to www.germanpool.com for the most up-to-date version

Product Structure



2. Control Panel

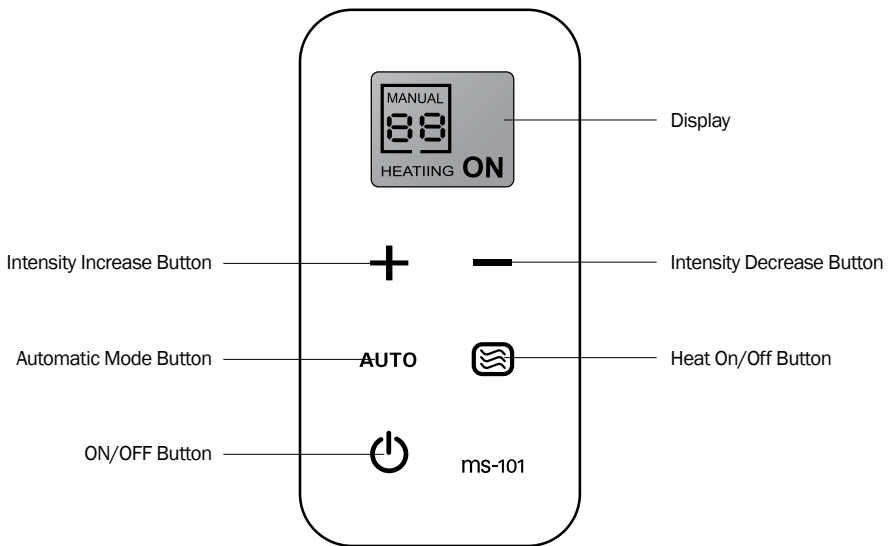*Product specifications are subject to change without prior notice. Please refer to www.germanpool.com for the most up-to-date version.

Operating Instructions

A) How to wear

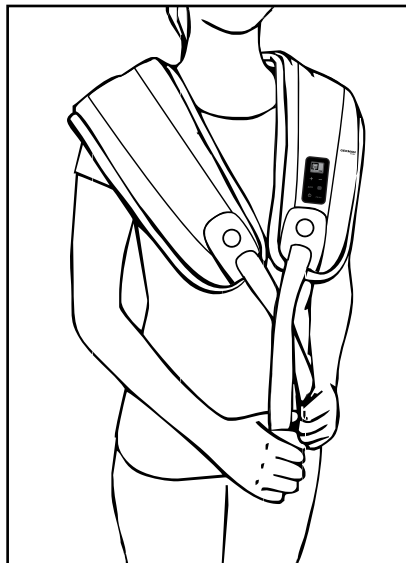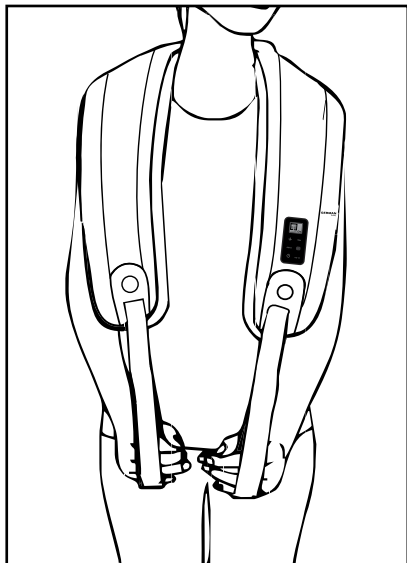With the control panel facing outward, apply the product over the body part to be massaged. You may hold the handle strap to secure the position or fasten the button on the handle straps.


Neck and shoulder


Back

Lower back

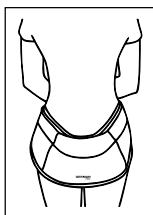Waist

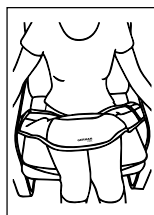Thighs

Operating Instructions

B) Operation

Note:

Make sure the unit's rated voltage matches the AC power before connecting the power plug. Place the massager on your body and make sure it is comfortably secured.

Manual Mode

1. Connect the appliance to a power socket. The unit will beep four times to indicate it is properly connected and ready for operation.
2. Press the ON/OFF Button  to start the unit. The unit will beep and begin the tapping massage. You can press the Intensity Increase/Decrease Buttons **+** , **-** to adjust the strength of the tapping. Each press of the Button is accompanied by a beep. There are twelve intensity levels in total, ranging from 1 to 12, and they are displayed correspondingly on the Display. If the Intensity Increase/Decrease is at a maximum, the unit will beep twice to indicate it cannot further increase/decrease. The start-up manual mode is preset at a 10-minute cycle with a tapping intensity of 1.
3. To turn on FIR heating, press the Heat On/Off Button . The appliance will beep once, and the Display will show "ON". The level of FIR heating cannot be adjusted; if the Button is pressed again, the unit will beep thrice and "ON" will disappear from the display to indicate that the FIR heating function is turned off.
4. To turn off the unit, press the ON/OFF Button . The appliance will beep four times and then turn off.

Automatic Mode

1. Connect the appliance to a power socket. The unit will beep four times to indicate it is properly connected and ready for operation.
2. Press the ON/OFF Button  to start the unit. Press the Automatic Mode Button . The unit will beep once and start in automatic mode; to exit automatic mode, press the icon again (the unit will beep three times).
3. The appliance will proceed to cycle through all 26 modes, automatically adjusting intensity from 1-12. The default time for this cycle is 10 minutes.
4. To turn on the FIR heating function, press the Heat On/Off Button . The appliance will beep once, and the display will show "ON". The level of FIR heating cannot be adjusted; if the Button is pressed again, the unit will beep thrice and "ON" will disappear from the display to indicate that the FIR heating function is turned off.
5. To turn off the unit in advance of the standard 10-minute cycle, press the ON/OFF Button . The appliance will beep four times then turn off.

Care & Maintenance

After use, and before cleaning this product, turn off the massager and unplug the power cord via the plug - not the cord. Pulling on the power cord may cause damage to the unit.

Do not use this product when in the bath or under water. The appliance is not meant for use in a damp environment.

If this product falls into water by accident, please turn off the power and unplug the power cord immediately. DO NOT remove the product directly from water.

Clean the unit with a dry cloth. DO NOT clean with a damp cloth, or with benzene-, thinner-, or any other solvent-based cleaning products on your massager and its parts.

Do not store the appliance near heat or open flame. Do not leave the massager in direct sunlight for a long period of time.

Technical Specifications

Model No	MS101
Voltage / Frequency	220 V ~ 50 Hz
Rated Power	80 W

* Product specifications are subject to change without prior notice. Please refer to www.germanpool.com for the most up-to-date version

* If there is any inconsistency or ambiguity between the English version and the Chinese version, the Chinese version shall prevail

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電郵：repairs@germanpool.com

Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online: www.germanpool.com/warranty.
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
 - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
 - damages caused by accidents of any kind (including transport and others);
 - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
 - purchase invoice is modified by unauthorized party;
 - product is used for any commercial or industrial applications;
 - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
 - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

Customer Service & Repair Centre

German Pool (Hong Kong) Limited

8/F, China Travel (Cargo) Logistics Centre
1 Cheong Tung Road

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This warranty card is valid only in Hong Kong

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© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

www.germanpool.com